

MVIEWM.005A



71308
10/13/01
PATENT

IN THE UNITED STATES PATENT AND TRADEMARK OFFICE

Applicant	:	Sherman, et al.	Group Art Unit 1652
Appl. No.	:	09/501,730	
Filed	:	February 10, 2000	
For	:	AGGREGATE-FREE URATE OXIDASE FOR PREPARATION OF NON- IMMUNOGENIC POLYMER CONJUGATES	
Examiner	:	Pak, Y.	

DECLARATION OF L. DAVID WILLIAMS, Ph.D.
SUBMITTED UNDER 37 C.F.R. § 1.131

Assistant Commissioner for Patents
Washington, D.C. 20231

Dear Sir:

1. I, L. David Williams, Ph.D., am a named inventor on the above-captioned patent application.
2. I have read and understood the specification and the claims of the above-captioned patent application.
3. I have read and understood Caliceti et al. *Bioconjugate Chem.* 1999, 10, 638-646, which has a web publication date of June 2, 1999.
4. Exhibit A, attached hereto, is a copy of several pages of one of my laboratory notebooks.
5. Pages 1 and 2 of Exhibit A show the results of a purification procedure that I conducted prior to June 2, 1999. The tracings shown on pages 1 and 2 are ion-exchange and size-exclusion chromatograms, respectively, of the product of my purification steps. The areas

Appl. No. : 09/501,730
Filed : February 10, 2000

under the peaks of the chromatograms on pages 1 and 2 show the relative abundance of several aggregates of uricase.

6. The areas under the peaks assigned to tetramers are much larger than the area under the peaks assigned to octamers or other aggregates. These data demonstrate the purification of uricase that is substantially free of aggregates larger than tetramers.
7. Therefore, prior to June 2, 1999 I had invented and was in possession the subject matter claimed in the above-captioned patent application.
8. I declare that all statements made herein are true, and that all statements made upon information and belief are believed to be true, and further, that these statements were made with the knowledge that willful, false statements and the like so made are punishable by fine or imprisonment, or both, under 18 U.S.C. § 1001, and that willful, false statements may jeopardize the validity of the application, or any patent issuing thereon.

October 4, 2001

Date

L. David Williams

L. David Williams, Ph.D.

S:\DOCS\SKT\SKT-1775.DOC
100301

linear5_M1_UV_280nm_01 — linear5_ConcB_01 — fractions_1

19-52

AU

Edit - Run Description

METHOD

Method base : ml

***** Main method *****

0.00 CONC_B	0.0	()
0.00 FLOW	1.0	()
0.00 RECORDER_SPEED	0.50	()
130.00 CONC_B	0.0	()
250.00 CONC_B	100.0	()
270.00 CONC_B	100.0	()

Sample BTG-Ply-KS-dt, 80%, 0.2u filtered

Loop load thru pump

Column MuxQ H210/10

Eluent A 0.7M pH 10.3 NaCO3

Eluent B 0.6M NaCl in Buffer A

Gradient Linear 0-100% C from 130 to 250 ml

Flow rate 1 mL/min

Fraction size 4 mL

RUN LOG

Method file name LINEARS.FMT

Run by LDW 8:52:26 am

Config: Method base is ml

Config: Injection valve is valve No. 1

Config: Pump A/B is P-500, calibration 110.00

Config: Pump C is P-1, calibration 150.00

Config: UV-M -> Monitor 1, Input range = 100 mV

Method start. Result name LINEARS.FRS

0.00	Auto:	Tube No 1	
0.00	Auto:	Pause	10:13:02 am
80.50	Auto:	Continue	10:16:54 am
130.02	Auto:	Tube No 2	
134.02	Auto:	Tube No 3	
138.01	Auto:	Tube No 4	
142.02	Auto:	Tube No 5	
146.01	Auto:	Tube No 6	
150.01	Auto:	Tube No 7	
154.01	Auto:	Tube No 8	
157.98	Auto:	Tube No 9	
161.99	Auto:	Tube No 10	
165.99	Auto:	Tube No 11	
169.97	Auto:	Tube No 12	
170.05	Manual:	Hold	11:46:34 am
173.96	Auto:	Tube No 13	
177.96	Auto:	Tube No 14	
181.95	Auto:	Tube No 15	
182.03	Manual:	Continue	11:58:33 am
185.96	Auto:	Tube No 16	
189.96	Auto:	Tube No 17	
191.97	Manual:	Hold	12:08:30 pm
193.95	Auto:	Tube No 18	
197.94	Auto:	Tube No 19	
201.94	Auto:	Tube No 20	
205.94	Auto:	Tube No 21	
208.00	Manual:	Continue	12:24:32 pm
209.94	Auto:	Tube No 22	
213.93	Auto:	Tube No 23	
217.92	Auto:	Tube No 24	
221.93	Auto:	Tube No 25	
225.93	Auto:	Tube No 26	
229.92	Auto:	Tube No 27	
233.91	Auto:	Tube No 28	
237.91	Auto:	Tube No 29	
241.90	Auto:	Tube No 30	
245.90	Auto:	Tube No 31	
249.89	Auto:	Tube No 32	
253.89	Auto:	Tube No 33	
257.87	Auto:	Tube No 34	
261.86	Auto:	Tube No 35	
265.87	Auto:	Tube No 36	
269.87	Auto:	Tube No 37	
273.85	Auto:	Tube No 38	
277.85	Auto:	Tube No 39	
281.84	Auto:	Tube No 40	
285.84	Auto:	Tube No 41	
298.02	Auto:	Method end	1:54:38 pm

1.5

1.0

0.5

0.0

0

50

100

12

13

14

15

16

17

18

19

20

21

22

23

24

25

26

27

28

29

30

31

32

33

34

35

36

37

38

39

40

41

42

43

44

45

46

47

48

49

50

51

52

53

54

55

56

57

58

59

60

61

62

63

64

65

66

67

68

69

70

tetramers_{us}octamer_{us}aggregates_{us}pool_{us}

Continued on Page 53

EXHIBIT A PAGE 1 OF 4

Read and Understood By

Mark Sifer

Signed

Date

Signed

Date

Flush Mono Q HR10/10 column in FPLC system with 0.6 M NaOH in 9 mM NaCl, then water, then a blank run.

Buffer A = 0.1 M pH 10.3 sodium carbonate, Buffer B = 0.6 M NaCl / buffer A.

Start loading 0.2 μ filtered BTG-Pig-KS-dN, lot AO-58-034 rec'd ~~from water-washed Nalgene~~ gross filtered wt 159.6 g
filter unit, load 90 mL @ 1 mL/min net wet wt 62.4 g

directly then Pump A. About 8.5 mL remaining

8:39:50 start load

8:44:14 pause at end of 1st piston cycle.

8:52:04 start Method LINEAR 5

10:12:38 pause method, stop load, start Buffer A


10:16:30 @ 80.5 mL continue

10:56:05 @ 120 mL

11:06:08 Start fraction collector, 4 min/fraction = 4 mL/fraction, at #2.
Hold gradient for three fractions at ~~41.6%~~ 33.3% B = 0.2 M NaCl
Hold gradient for 4 fractions at 41.6% B = 0.25 M NaCl

13:34:14 reach 100% B, 250 mL

13:54:14 pump off, end @ 270 mL

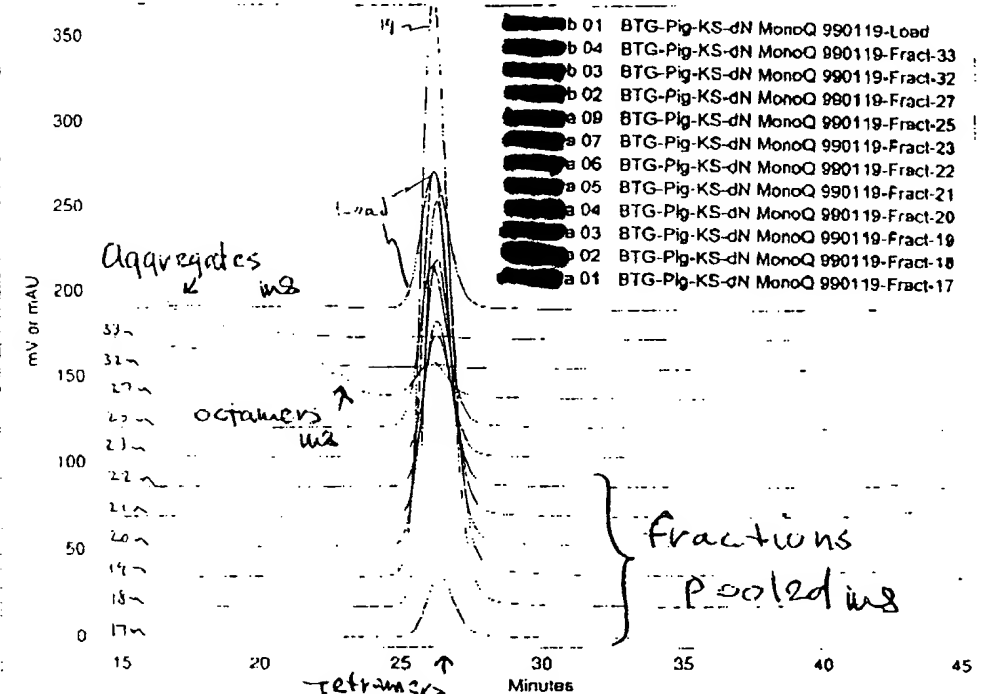

Superdex 200 HR10/30
Buffer: 10 mM Carbonate + 100 mM NaCl, pH 10.1; 0.5 mL/min
UV2000 B 214nm = UV276nm

Continued on Page 54

EXHIBIT A PAGE 2 OF 4

Read and Understood By

Mark Sifer

Signed

Date

Signed

Date

Remove 180 μ l returns of 19-53-Fract's 17-22 into HPLC vials w inserts. Also prepare dilutions as follows

17	1x	360 μ l return
18	1/8	50 μ l + 350 μ l buffer ^{water} low
19	1/12	30 μ l + 330 μ l water
20	1/8	50 μ l + 350 μ l water
21	1/5	70 μ l + 280 μ l water
22	1/3	100 μ l + 200 μ l water

Pool remaining 17-22 after return using fract 17 as wash of other to Pool wt = 21.565g

Remove 2 mL Pool return
 $(21.565 - 2) \times 10.6 \frac{\text{mg}}{\text{mL}} = 207 \text{ mg target}$

Naigh out 233 mg NPC-10K-PEG (Skanwater LG-128-01)

and add 0.932 mL 1 mM HCE at 18:16. Momentary yellow color on solids before they are well suspended. residual base in NPC-PEG?

Vt of added HCE = 0.941g. Add 1.099g PEG soln dropwise with stirring to the 19.5 mL Pool, ice bath, over 1-2 min. 1 mL PEG soln weighed 1.028g. Calcd NPC PEG added = $233 \text{ mg} \times \frac{1.099}{0.941 + 233} = 218 \text{ mg}$
 Covered + refrigerated

Remove μ l sample of rxn for HPLC

Prepared 10K Midbee ultrafilter

by rinsing with water, then 0.1 M NaOH

Rinsed with water, then 0.1 M pH 10.3 buffer (left in lines, i.e. ~3 mL) and began ultrafiltering round 1 reaction at 16:10. By 16:23 of about 5 permeate had accumulated. Continue concentrating to minimum vol. Stop at 17:37. Remove 3.266g (~3.15 mL) retentate. Rinse apparatus with 6 mL water, save separately.

Begin stirring retentate in ice bath. Remove 110 μ l sample + dilute with 140 μ l (1/20) carbonate HPLC buffer at ~10:15. NPC-10K PEG weighed yesterday (344 mg initially, dropping to 337 mg over 1-2 min static electricity?) + stored covered in freezer. Added 337 mg x 4 μ l x 90% = 1.22 mL 1 mM HCE at 10:15. Again, yellow color on solid until well suspended. Refrigerate until add to rxn dropwise (~1 drop/sec) 10:54:23 - 11:04:10. Rinse tube with 0.12 mL mM HCE after reducing foam under partial vacuum. Add rxn to rxn at 11:07. Add 2 mL carbonate HPLC buffer to PEG-tube. $UV_{400 \text{ nm}} @ 13:38 = \frac{1.01}{18.4 \text{ g} \cdot \text{cm}^{-1} \cdot \text{mm}} = 0.055 \text{ mm} \times 10 \frac{\text{mg}}{\text{mL}} = 0.55 \text{ mg}$
 At 11:18:00 remove 20 μ l rxn sample for W monitoring at room temp. $t_{1/2} = 11.5 \text{ min}$ (see LDW 01214.KD). Drop 0.1 mm W cell, breaking it. Dilute ~~10 μ l~~ UV soln all recovered (<20 μ l) UV soln with 380 μ l water (~1/20) for HPLC 9901214.01. PEG submit calcd as 7.5, so

Continued on Page 5

Read and Understood By

L D Williams

Signed

Date

Mark Sifer

Signed

Date

PROJECT _____

quench 2 mL in 22 mL 0.60 M Gly Gly ()
Filter through 0.45 μ filter. See [redacted] & OI for HPLC.

~~Start ultrafiltration on 30K Membrane~~

Took wrong tube from refrigerator: ultrafiltering Round 1 permeate
No activity in retentate or permeate. Ultrafiltered 0.1 mg/mL BSA/50 mM
brax to try to chase off bound activity before discovering there was no
activity to chase off.

[redacted] Round 3: Removed 10 μ L sample of unquenched round 2 rxn
for diln with 190 μ L water + HPLC.

Dissolved 102.4 mg NPC-10K-PEG in 0.36 mL 1 mM HCl and added
to the round 2 rxn at 17:23 with stirring in an ice bath. Rinsed in
PEG with 40 μ L 1 mM HCl.

EXHIBIT A PAGE 4 OF 4

Continued on Page

Read and Understood By

Mark Sawyer

Signed

Date

Signed

Date